

as per *Dr. Mihaljević*



HEAT-INSULATING MATTRESS FOR:

as per Dr. Mihaljević


- DENTAL TREATMENTS
- DENTAL X-RAYS
- SURGICAL PROCEDURES
- X-RAYS
- CT, MRT
- SONOGRAPHY
- CAGES
- TREATMENT TABLES

DENTO-DORM SITTING, LYING AND STORAGE MATTRESSES

Developed, produced and clinically tested in Germany and Austria, according to strict standards. Due to the flexible positioning wedges as all as the anatomically shaped head-positioning wedges, the comfort of the veterinarian and the patient are improved in the surgery and X-ray area. The correct positioning of the patient increases the quality of the results, for example of dental X-rays.

POSITIONING MATTRESS FOR OPERATIONS, DENTAL TREATMENTS, X-RAYS AND SURGERY

Hygiene due to a sealed surface

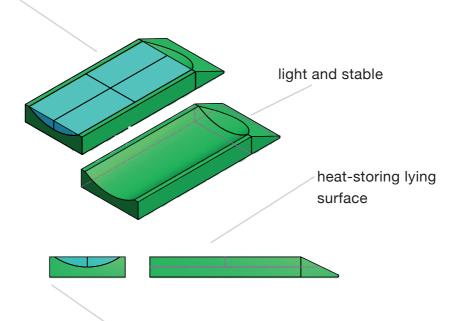

no shifting on slippery surfaces and diagonally set tables, adjustable through 4 positioning wedges



positioning mattress "standard" (black)



positioning mattress CT/MRT (silver grey)



FEATURES OF THE DENTO-DORM® POSITIONING MATTRESSES

DENTO-DORM®

• stabilizes the body temperature of the patient during anaesthesia thanks to numerous air cushions.

DENTO-DORM®

• prevents the patients from shifting even when in a sloping position.

The positioning mattress M "standard" includes 4 additional positioning wedges.

DENTO-DORM®

• supports the positioning of the patients, thereby enabling ergonomic working conditions.

DENTO-DORM®

• available in various sizes suitable to position any type of pet including large dogs.

SAFE POSITIONING AND HEAT-INSULATION OF THE PATIENTS!

DENTO-DORM®

• water-proof thanks to a hygienic, small pore, sealed surface which is easy to clean and disinfect.

DENTO-DORM®

 does not contain any chemical solutions, is robust, stable and very light.

DENTO-DORM®

• has a concavely swung shape which enables a comfortable and secure positioning of the patients and provides a special positioning aid during dental X-rays.

DENTO-DORM®

• various uses are possible thanks to the two differently designed ends.



INNOVATIVE PRODUCTS FOR THE VETERINARIAN FIELD



DENTO-DORM®

TREATMENT TABLE MATTRESS

Robust polyethylene core. Anti-skid, -pierce, -scratch and solvent resistant cover, welded to be water-proof. Easy to clean and disinfect. Colour black, red, green. Clinically tested. (Tested by LABOKLIN).

Standard size: fits tables

Size 1 1300x500 mm

Art-No. 29767-032

Size 2 1300x600 mm

Art-No. 29767-031

Other sizes of all table manufacturers available.

*Note: anti-parasitic spot-on preparations can lead to discolouration; this does not pose any kind of impairment of the material's features.



POSITIONING MATTRESS Following proper disinfection, bacteriologically and virologically safe.

LABOKLIN

DENTO-DORM®

POSITIONING MATTRESS "standard"

Made from hard foam, hygienic, small pore, sealed surface. Robust and light material. Waterproof, can be wiped and disinfected. Clinically tested. Bacteriologically and virologically harmless after proper disinfection. (Tested by **LABOKLIN**). Incl. 4 positioning wedges.

Area of use: Dental treatment, X-rays, sonography, blood collection.

Size M 750x300x80 mm

Art-No. 29767-001 / suitable for cats and dogs

Size L 1000x350x120 mm

Art-No. 29767-003 / suitable for large dogs





DENTO-DORM®

HEAD POSITIONING WEDGES - SET

Made from hard foam, hygienic, small pore, sealed surface. Robust and light material. Waterproof, can be wiped and disinfected. Clinically tested. Bacteriologically and virologically harmless after proper disinfection. (Tested by LABOKLN). The set includes 2 head positioning wedges which fit the DENTO-DORM® positioning mattress size M and 1 head positioning wedge which fit the DENTO-DORM® mattress size L.

Set M -2 pieces-

Art-Nr. 29767-006

Set L -1 piece-

Art-Nr. 29767-007



DENTO-DORM POSITIONING MATTRESS Following proper disinfection, bacteriologically and virologically safe.

DENTO-DORM®

POSITIONING MATTRESS "standard"

Made from hard foam, hygienic, small pore, sealed surface. Robust and light material. Waterproof, can be wiped and disinfected. Clinically tested. Bacteriologically and virologically harmless after proper disinfection. (Tested by LABOKLIN). Incl. 2 positioning wedges as well as 1 head positioning wedge with 80 mm and 1 head positioning wedge with 120 mm. Area of use: Dental treatment, X-rays, sonography, blood collection.

Size S 480x200x80 mm

Art-No. 29767-002 / suitable for pets

INNOVATIVE PRODUCTS FOR THE VETERINARIAN FIELD



DENTO-DORM®

CAGE MATTRESS

2-layered mattress core made from hard and soft foam with a PVC cover, as well as 100% waterproof seams. Two-coloured light grey/green. Easy-care, scratch-free mattress which can be wiped and disinfected with two functionally different sides:

- the light grey side is pleasantly soft (optimal sleeping comfort for lying patients)
- the green side is made of a firmer material (for patients who move about the cage), has the positive effect of a cosy, heat-storing sleeping mattress.

Fits in Shor-Line cages.

Size 1 700 x 540 x 60 mm

Art-No. 29767-050

Size 2 1140 x 680 x 60 mm

Art-No. 29767-051

Other sizes available on request.



DENTO-DORM®

POSITIONING MATTRESS CT/MRT

Made from hard foam, hygienic, small pore, sealed surface. Robust and light material. Colour silver grey. Waterproof, can be wiped and disinfected. Dermatologically safe, clinically tested. Incl. 4 positioning wedges, 1 head positioning wedge and 2 special stabilizing wedges for the limbs. The special shape of the super light and heat-storing DENTO-DORM® positioning mattress CT/MRT and the additional positioning wedge enables a secure and stable positioning of the patient during CT, MRT, ultrasound or X-ray examinations.

Size S (750 x 270 x 100 mm)

Art-No. 29767-047 / suitable for pets

Size M (1000 x 350 x 150 mm)

Art-No. 29767-048 | suitable for dogs and cats



DENTO-DORM®

TRACHEAL-SPONGE

Made from cellular rubber, this is a *sealing sponge* for the trachea, especially during calculus removal and dental treatments. The DENTO-DORM® Tracheal-Sponge has a tightly sealing outlet opening for the tracheal tube. Available as a basic case kit with various Tracheal-Sponge sizes for small cats to large dogs. The case comes with 48 Tracheal-Sponges. Bacteriologically and virologically harmless after proper disinfection. (Tested by LABOKLIN).

Tracheal-Sponge Case Art-No. 29767-045 48 pieces Tracheal-Sponges.

Refill -rough- Art-No. 29767-0451 12 pieces Tracheal-Sponges.

Refill -fine- Art-No. 29767-0452 36 pieces Tracheal-Sponges.

EXAMPLES OF USE

DENTO-DORM®

POSITIONING MATTRESS "Standard"



Positioning mattress Area of use: dental treatment, X-rays, sonography, blood collection





Waterproof, light, body-heat-storing mattress. Does not slip from the sloped table. Simplifies standardised X-raying by precisely positioning the patient. Also makes ergonomic operating easier, while at the same time offering the patient more comfort.

DENTO-DORM®

HEAD POSITIONING WEDGE



Head Positioning Wedge Set





Enables an exact positioning of the head for dental X-rays, CT and surgeries. Recesses for the limbs on the bottom. Fits on the positioning mattress size M or size L.

DENTO-DORM®

POSITIONING MATTRESS CT/MRT



Area of use: CT/MRT, X-rays, sonography, surgery



Prone position



Back position



CT/MRT



CT/MRT Area

Thanks to its concave shape, this mattress is excellent to position patients, including narrow-chested animals, for X-rays, CT, MRT, surgeries, sonographies.

EXAMPLES OF USE

DENTO-DORM®

TREATMENT TABLE MATTRESS



Treatment table mattress

Robust, long-lasting table-top. Mattress made of scratch- and pierce-free material which does not slip away on smooth surfaces - even for heavy, restless dogs. Clinically tested. A pleasant table-top for the animals, which results in patients being less stressed during examinations and treatments.



New feature: 100% waterproof welded without seam

DENTO-DORM®

CAGE MATTRESS



The green side is made of a firmer material (for patients who move around the cage)



Cage mattresses with the unique feature of two different functional sides, which are indicated by their colour. They are gentle on the joints, especially for older and heavier patients.

All sizes can be requested.



The light grey side is pleasantly soft (optimal sleeping comfort for lying patients)



Optimal sleeping comfort for lying patients after surgery, raises temperature and allows for a pleasant stay for patients in a cage.

EXAMPLES OF USE

DENTO-DORM®

TRACHEAL-SPONGE





Dental treatment



Sponge seals with central, densely sealing perforations for the tracheal tube. Employed during calculus removal and dental treatments far in the pharynx and prevents calculus fragments and water from reaching the pharynx. At the same time ensures the tube from slipping.

The material resists commercial disinfectants and autoclaving, so that - after used properly - multiple use is bacteriologically and virologically harmless.

Tested by **LABOKLIN**. Detailed cleaning tips are included in every case.



Application tips can be found on our YouTube-Channel: Dento-Dorm.

OPERATING AT A HIGHER LEVEL

The idea was born during the daily routine in surgery. On the one hand, it was meant to offer the animal optimal heat-insulating lying and sleep comfort. On the other, the treating veterinarian's work was supposed to be made easier and more ergonomic, improving the results at the same time. This resulted in the completely newly developed, high quality and long-lasting DENTO-DORM® products.

PROTECTED SYSTEM

The system which was designed by us was granted the utility model protection by the German Patent and Trade Mark Office on February 10, 2011 (No. 20 2010 014 780.8).

SECURED QUALITY

DENTO-DORM® is developed, produced and clinically tested in Germany and Austria according to strict standards. Thanks to physically interlaced and stable foam material with a sealed surface, we can offer our customers the highest quality.

SOPHISTICATED DESIGN

Our heat-insulating positioning mattresses offer highest flexibility for dental treatments, X-rays, sonographies, blood collection and surgeries on dogs, cats and small pets. They can be combined in all variations.

INNOVATIVE PRODUCTS FOR THE VETERINARY SECTOR & APPLICATION EXAMPLES

DENTO-DORM®

LUXATOR RACIFEIN as per Dr. Mihaljević

The 'original' with absolutely straight, delicate and sharp handling tip to luxate small teeth with a narrow periodontal gap and to remove root splinters and small bone fragments. With the help of these finely crafted dental instruments, luxating and extracting especially cat teeth and small dog teeth becomes much easier, as the very sharp and thin handling tip allows access even to narrow periodontal gaps.

To do so, the Luxator is inserted with apical pressure around the tooth and in the imaginary line of the tooth root (i.e. very steeply and not at a 45-degree angle to the tooth). Following that, slightly rotating movements from the wrist are made as well as a further deeper insertion of the Luxator with apical pressure and subsequent slight rotation. This procedure is performed several times and always on different sides of the tooth until the tooth is mobilised.

The sharp tip works like a cutting tool to cut through the circular periodontal fibres in the gap. Regular sharpening is thus recommended. While handling the tool, always make sure that the instrument is guided securely (fingers far down towards the handling tip) and provide proper support to prevent stabbing injuries due to accidental slipping.







DENTO-DORM®

LUXATOR RACIFEIN 1 MM as per Dr. Mihaljević

Very gracile handling tip, particularly suitable for luxating small cat teeth (Incisive, P2 or broken-off root remains)

Handling tip: 1 mm wide, 10 mm long

Total length: 14 cm

Size 1 (1 mm w x 10 mm l) Total length 14 cm

Item no. 29767-100



DENTO-DORM®

LUXATOR RACIFEIN 2 MM as per Dr. Mihaljević

Delicately manufactured handling tip, excellently suited for luxating cat and small dog teeth.

Handling tip: 2 mm wide, 12 mm long

Total length: 14 cm

Size 2 (2 mm w x 12 mm l) Total length 14 cm

Item no. 29767-101



DENTO-DORM®

Luxator Racifein 1 mm as per Dr. Mihaljević





DENTO-DORM®

Luxator Racifein 2 mm as per Dr. Mihaljević





as per Dr. Mihaljević



DENTO-DORM GbR

Im Andermannsberg 46 88212 Ravensburg

www.dento-dorm.de info@dento-dorm.de





